

# One Standard Across the Care Continuum

BSM-6000 series  
bedside monitors



# Breaking Down Barriers to Care

The Institute of Healthcare Improvement's triple aim—designed to push quality initiatives while simultaneously fostering marketplace efficiencies—has rippled across the healthcare continuum. Advancement simply for the sake of change is not enough; it must enhance healthcare in some way. Improvements must be quantifiable, as well as qualifiable. Care must be coordinated. Every clinician must manage greater risk as patients present with multiple comorbidities. To navigate this dynamic, providers across the spectrum—from the smallest practices to the largest health systems—need allies capable of helping them remain ahead of the rapidly shifting environment while ensuring the highest quality of care for patients.



A soft-focus photograph of a medical setting. On the left, a woman in a white lab coat with a stethoscope around her neck is looking towards the right. On the right, a patient is lying in a hospital bed, looking up at the doctor. The patient is wearing a teal hospital gown. The background is bright and out of focus, showing what appears to be a window with sheer curtains.

# Nihon Kohden is that ally.

A leader in precision medical products and services, only Nihon Kohden offers reliable, integrated multi-modality products that serve patients across all care areas. We bring quality clinical solutions that provide access to a deeper, more comprehensive level of information, enabling more accurate diagnoses and ultimately, better outcomes.

# Premium-as-Standard Design

We believe every patient deserves the highest standard of care. Nihon Kohden's premium-as-standard philosophy is the belief that every monitor should be fully appointed with features—both standard and premium—unlocked and ready to use at a moment's notice. This ensures that our technology can be employed in the broadest range of acuity levels and seamlessly transition between care areas as patient need dictates.

The BSM-6000 series of monitors offer unrivaled technology that works across the healthcare continuum and realizes our premium-as-standard philosophy, allowing providers to deliver care without compromise.

Some of our premium-as-standard capabilities include:

- Comprehensive arrhythmia detection and recall, including advanced Atrial Fibrillation algorithm
- Multi-waveform/multi-parameter full disclosure
- ST segment analysis as well as diagnostic 12 lead ECG capability
- Drug, hemodynamic and pulmonary calculations



**BSM-6301**  
10.4-inch LCD



**BSM-6501**  
12.1-inch LCD

Nihon Kohden's unique Smart Cable™ technology miniaturizes circuits found in traditional modules and embeds that circuitry into a Smart Patient Cable. When you plug a Smart Cable into a Multiport, the associated parameter is automatically detected, displayed and measured. With this technology, you'll get parameter flexibility at a significantly reduced cost with seamless and immediate access to blood pressure, cardiac output, EtCO<sub>2</sub>, temperature, BIS and more, when and where you need it for rapid clinical assessment across care areas.




01  
CD

BSM-6701  
15-inch LCD




## Arrhythmia

The BSM-6000 Series provides high accuracy multi-lead arrhythmia detection and storage of over 16,000 arrhythmia events. Each event is time-linked to the full disclosure waveforms to determine what led up to, and what followed, the captured event.

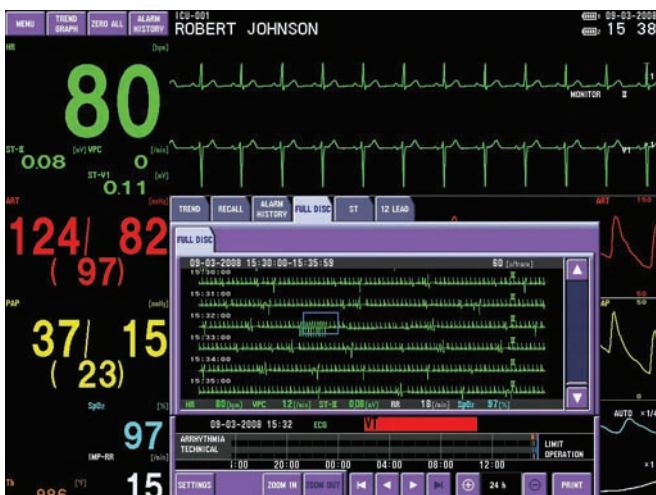

## Full Disclosure

Full disclosure waveforms allow you to validate alarm and numeric findings and to make treatment decisions based on more accurate monitored data. The BSM-6000 Series provides storage and review capabilities at the bedside monitor that are typically found only in a central station.

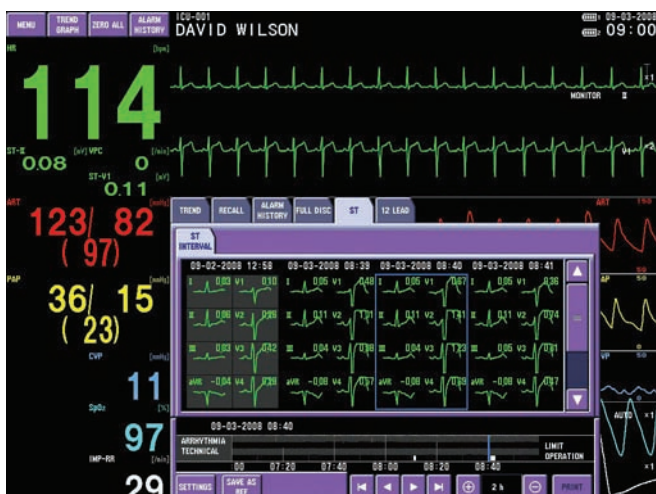

## ST Template

Multi-lead ST segment monitoring provides you with continuous oversight to transient changes in your patients' cardiac condition and are stored minute-to-minute in the monitor for comparison.

# Redefining Transport

Continuity of monitored care during transport ensures the highest quality of care for patients. The unique design of Nihon Kohden's BSM-1700 Transport monitor and the BSM-6000 redefines transport monitoring. Simply disconnect the BSM-1700 Transport monitor from the BSM-6000 monitor or Data Acquisition Unit and your patient can be transported with all monitoring capabilities remaining the same. When the patient is transferred to their new care setting and the BSM-1700 is reconnected to another BSM-6000 or Life Scope G9 monitor and patient information, including full disclosure, is uploaded to the new bedside display creating one seamless, reliable patient record enhancing workflow and care coordination.



One additional benefit of the Data Acquisition Unit is that it can be extended via an umbilical cord to be located on a bed rail, gurney or IV pole next to the patient. This alleviates the hassles associated with cable management at the bedside. Since the Data Acquisition Unit contains user function keys, it can be placed on either side of the patient for optimal clinical workflow efficiency.



 **NIHON KOHDEN**  
[us.nihonkohden.com](http://us.nihonkohden.com)

MLLB 011 (B)-CO-0511

Different Thinking for  
Better Healthcare.™

Different Thinking for Better Healthcare is a trademark of Nihon Kohden.