

Data at Your Fingertips

- Provides a secure method for viewing a range of Nihon Kohden-generated patient data that empowers clinicians to make timely decisions and supports quality of care
- Offers a continuous and simultaneous display of data for up to 12 patients on iPad[®] devices and up to eight patients on iPhone[®], Google Pixel[™] 1, Samsung Galaxy[®] Note8, and Zebra[®] TC51 Android[™] devices
- Offers a robust, easy-to-use interface, which allows users to see current waveforms, vital signs, stored data, and much more
- Compatible with Apple mobile iOS devices and Android devices

Interactive Monitoring Data

- Graphical and tabular trends
- Hemodynamic calculations
- Arrhythmia events
- Minute-to-minute ST templates
- Multi-parameter full disclosure waveforms
- Interpretative 12-lead ECGs


Secure, Real-Time Viewing

- Secured access for local or remote users
- Near real-time waveforms and numerics with user-selectable settings for customized views
- Allows data to be viewed in near real-time within the hospital network or remotely via a VPN connection
- Leverages NK-HiQ[™] Enterprise Gateway multi-level security controls, which restrict access to patient data based on user, user groups, and/or mobile client



ViTrac

Remote Viewing App

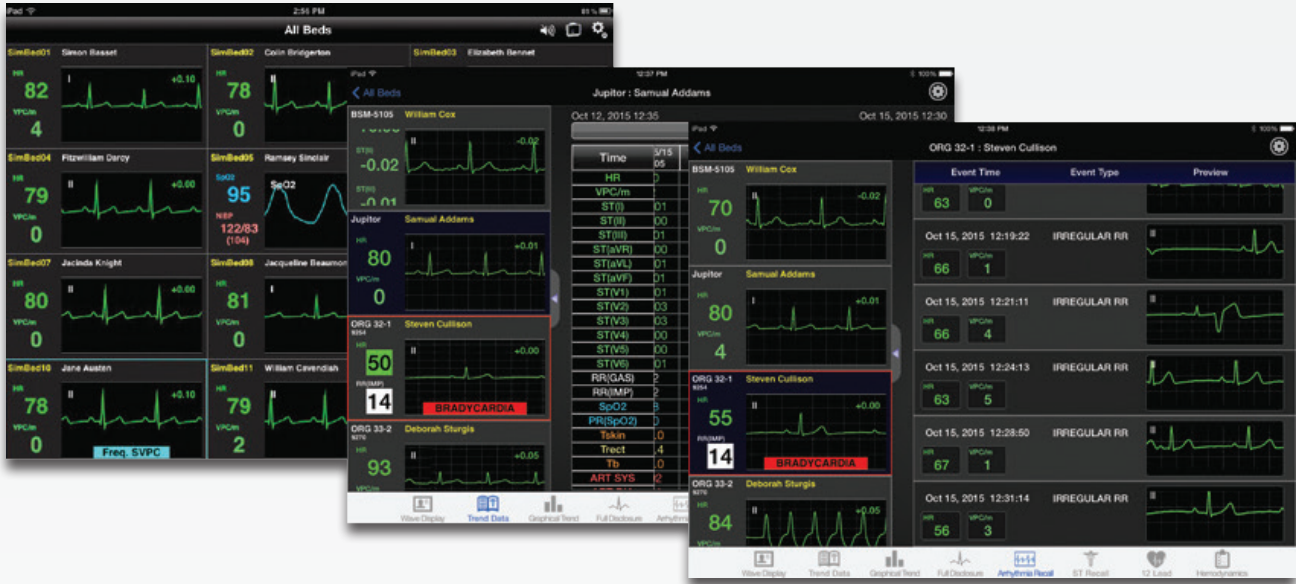

Compatible with Multiple Devices

- Apple iPad, 2nd generation or higher, running iOS 6.1 or higher
- Apple iPad mini, 1st generation or higher, running iOS 6.1 or higher
- Apple iPhone 4 or higher, running iOS 6.1 or higher
- Apple iPod touch®, 4th generation or higher, running iOS 6.1 or higher
- Google Pixel 1, Samsung Galaxy Note8, and Zebra TC51 running Android 4.1 (Jelly Bean) and higher, up-to Android 8.1 (Oreo)

System Overview

The ViTrac remote viewing application is part of the NK-HiQ Enterprise Gateway platform. Designed for secondary viewing, the application includes all required software to provide clinicians with secure remote access to patient data.

As the backbone of the Nihon Kohden Patient Monitoring System, the NK-HiQ Enterprise Gateway platform consolidates multiple applications into a single powerful server — saving time, space, and money.

For more information, please contact us at 1-800-325-0283 or visit us.nihonkohden.com

Different Thinking for Better Healthcare.®

